

Beds

– *Comfortable, practical and multi-adjustable*

Guldmann beds are built for maximum comfort and easy use, with a practical design that saves time, work and effort for all staff groups.

There are two series of Guldmann beds – one specifically for hospitals, and one for all other forms of care. Their unique design ensures effective care for the user, and enables the carer to work in the most ergonomically correct way. They are also easy for porters, technical staff and cleaners to move, service and clean effectively.

Guldmann™

Complete system

Guldmann beds are designed as an integral part of a well-thought-out infrastructure for all-round bed care, including lifting and moving people within nursing and care environments.

The entire bed is designed with this in mind, with a base divided into four adjustable sections. These can be raised and lowered in ways that provide the user with the best possible comfort and the most effective nursing care. For example, the thigh section can be set at angles of as much as 45 degrees to make sure the user does not slide down the bed, and has maximum leverage to help movement. In this position, the buttocks remain on the flat “comfort board”, and it is easier to deal with issues relating to pressure points.



Guldmann™

Guldmann beds serve as an active tool that helps ensure more time for care procedures. For example, a trained carer can combine different positions to help move the user around in the bed, using no physical effort at all. The sections can also be positioned to make it easy for the carer to undertake specific tasks, while the user remains comfortable.

The beds can be positioned to make it as easy as possible to place the user in a sling ready for lifting and moving, and transfer to other equipment and aids. Raising or lowering the bed also makes it easy for ambulant users to get in and out of bed with minimum assistance from the carer.



Multi-adjustable, practical

These beds are easy to adjust and move around, helping reduce stresses and strains by enabling the carer to work in the most appropriate position at all times. This makes these beds intensely practical within a professional care environment, where it is necessary to meet widely varying needs.

The remote control unit for Guldmann beds can be operated by both the user and the carer, and the lifting mechanism is smooth and virtually silent.

The beds are fully adjustable in height, and the four sections of the bed base can be combined into a vast range of positions. Electric motors are available to adjust each section independently, for maximum comfort and the most effective nursing care. As option there is battery back-up for use when the bed is not connected to mains power, for example while in transit.

Well-designed

Every design detail in these beds has been carefully considered to maximize the user's comfort and minimize the work effort of all staff. Guldmann bed frames are free of fittings and brackets, all operating handles are placed so that they are easy to get at, and adjusting or servicing these beds requires no tools of any kind.

If needed, bed guard rails are mounted on the frame itself, keeping the overall width of the bed to a minimum.

Guldmann beds can be placed on special trolleys for easy transport and space-saving storage when not in use.



Technical specifications

Guldmann™

SWL	Max user weight	Number of motors	Central brake system	Width cm		Item Nos.			
						Guldmann care bed	Guldmann hospital bed	Guldmann insertion bed	Guldmann bed lift
240 kg	200 kg	2		82	standard head/foot boards	P07-2U82SE			
240 kg	200 kg	3		82	standard head/foot boards	P07-3U82SE			
240 kg	200 kg	4		82	standard head/foot boards	P07-4U82SE			
240 kg	200 kg	2		90	standard head/foot boards	P07-2U90SE			
240 kg	200 kg	3		90	standard head/foot boards	P07-3U90SE			
240 kg	200 kg	4		90	standard head/foot boards	P07-4U90SE			
240 kg	200 kg	2	C	82	standard head/foot boards	P07-2C82SE			
240 kg	200 kg	3	C	82	standard head/foot boards	P07-3C82SE	H07-3C82SE		
240 kg	200 kg	4	C	82	standard head/foot boards	P07-4C82SE	H07-4C82SE		
240 kg	200 kg	2	C	90	standard head/foot boards	P07-2C90SE			
240 kg	200 kg	3	C	90	standard head/foot boards	P07-3C90SE			
240 kg	200 kg	4	C	90	standard head/foot boards	P07-4C90SE			
<hr/>									
240 kg	200 kg	2		82	wooden head/foot boards	P07-2U82LE			
240 kg	200 kg	3		82	wooden head/foot boards	P07-3U82LE			
240 kg	200 kg	4		82	wooden head/foot boards	P07-4U82LE			
240 kg	200 kg	2	C	82	wooden head/foot boards	P07-2C82LE			
240 kg	200 kg	3	C	82	wooden head/foot boards	P07-3C82LE			
240 kg	200 kg	4	C	82	wooden head/foot boards	P07-4C82LE			
<hr/>									
240 kg	200 kg	3		82	without bed frame		I07-3U82XE		
240 kg	200 kg	4		82	without bed frame		I07-4U82XE		
240 kg	200 kg	3	C	82	without bed frame		I07-3C82XE		
240 kg	200 kg	4	C	82	without bed frame		I07-4C82XE		
<hr/>									
240 kg	200 kg				without bed board				S07-UE
240 kg	200 kg		C		without bed board				S07-CE



Identification of item Nos.

P07-	3	C	82	SE
P07 Guldmann care bed H07 Guldmann hospital bed I07 Guldmann insertion bed S07 Guldmann bed lift	2 motors: lifting motor back section motor 3 motors: lifting motor back section motor leg section motor 4 motors: lifting motor back section motor leg section motor thigh section motor	C with central brake system U without central brake system	width in cm	SE standard head/foot boards prepared for steel bed guards SE standard for hospital bed LE wooden head/foot boards with wooden bed guards XE without bed frame E without bed board

Beds

Guldmann™



2

Care bed, wood

Care bed with 2-4 motors with fixed wooden head/foot boards and wooden bed guards.
..... P07-



3

Hospital bed

Hospital bed with 3-4 motors and Trendelenburg function incl. steel bed guards and central brake system.
..... H07-



4

Bed lift

Intended for use with an ordinary bed, when there is a need for electric height adjustment.
..... S07-



5

Insertion bed

To be placed in an ordinary bed with motors for operation of back, thigh, and leg sections as well as for height adjustment.
..... I07-




1

Care bed

Care bed with 2-4 motors with detachable wooden head/foot board. The bed is prepared for steel bed guards.
..... P07-

Accessories for bed

Guldmann™

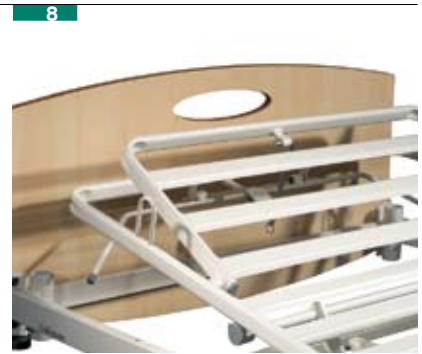

1
Bed guard, steel – standard
 right, care bed 4310-1
 left, care bed 4311-1
 right, hospital bed,
 chromium-plated, 4310-6
 left, hospital bed,
 chromium-plated, 4311-6

2
Bed guard, steel – with cord drive
 right, care bed 4310-2
 left, care bed 4311-2

3
Bed guard, steel – high
 right, care bed 4410-1
 left, care bed 4411-1
 right, hospital bed 4410-3
 left, hospital bed 4411-3

4
Bed guard, steel – small
 446000

5
Extension piece for bed guard, steel
 446100

6
Cover for steel bed guard
 43700

7
Cover for wooden bed guard
 43710

8
Head support
 for care bed 441600-1
 for hospital bed 441600-3

9
Support handle
 4478

Accessories for bed

Guldmann™



- 10**
Mattress
 82 x 200 x 13 cm, Multi flex 9158
 82 x 200 x 12 cm, cold foam 9183

- 11**
Mattress extension
 82 x 20 x 13 cm,
 for Multi flex mattress 9157
 82 x 20 x 12 cm,
 for cold foam mattress 9189

- 12**
Mattress cover
 Cotton 9187
 Incontinence 9188

- 13**
Lifting pole
 for care bed, max 90 kg 4320-1
 for hospital bed, chromium-plated,
 max 90 kg 4320-6
 for care bed, max 75 kg 4325-1

- 14**
Fender wheel
 horizontal 6178
 vertical 446200

- 15**
Bed trolley
 1050

- 16**
Bedside table
 right/left 685